

Sylwia BABICZ<sup>1</sup>, Artur ZIELIŃSKI<sup>2</sup>, Janusz SMULKO<sup>1</sup>, Kazimierz DAROWICKI<sup>2</sup>

<sup>1</sup>GDANSK UNIVERSITY OF TECHNOLOGY, FACULTY OF ELECTRONICS, TELECOMMUNICATIONS AND INFORMATICS

<sup>2</sup>GDANSK UNIVERSITY OF TECHNOLOGY, CHEMICAL FACULTY, 11/12 Gabriela Narutowicza Street, 80-233 Gdansk, Poland

## A measurement system for nonlinear surface spectroscopy with an atomic force microscope during corrosion process monitoring

M.Sc. Sylwia BABICZ

Ph. D. student at the Department of Metrology and Optoelectronics, Gdansk University of Technology, Poland. Her current research interests are focused on higher harmonic imaging by scanning tip of atomic force microscope.



e-mail: sylwia.babicz@eti.pg.gda.pl

Ph.D., D.Sc. Janusz SMULKO

Professor at the Department of Metrology and Optoelectronics Gdansk University of Technology, Poland. He is a member of the board of the Polish Society of Theoretical and Applied Electrotechnics in Gdansk, a member of the Microsystems and Sensors Section of the Metrology Committee of the Polish Academy of Science and a member of IEEE. Editor-In-Chief of Metrology and Measurement Systems Journal.



e-mail: jsmulko@eti.pg.gda.pl

Ph.D. Artur ZIELIŃSKI

Currently a research staff member and lecturer at the Chemical Faculty of Gdańsk University of Technology, where he received his PhD grade in chemical sciences. His current research interests are in the field of scanning probe microscopy and nanoimpedance spectroscopy.



e-mail: ziela@chem.pg.gda.pl

Prof. Kazimierz DAROWICKI

Titular professor since 2001, head of Department of Electrochemistry, Corrosion and Materials Engineering, Chemical Faculty of Gdańsk University of Technology since 1996. Area of scientific interests: techniques for characterization of electrodes and electrochemical processes, spectral analysis, corrosion and anticorrosion technologies. Author or co-author of over 300 journal papers, 6 textbooks and over 500 studies for industry.



e-mail: zak@chem.pg.gda.pl

### Abstract

In addition to traditional imaging of the surface, atomic force microscopy (AFM) enables wide variety of additional measurements. One of them is higher harmonic imaging. In the tapping mode the nonlinear contact between a tip and specimen results in higher frequency vibrations. More information available from the higher harmonics analysis proves to be helpful for more detailed imaging. Such visualization is especially useful for heterogeneous surfaces which are studied to understand corrosion mechanisms. In this paper the measurement system for nonlinear surface spectroscopy by AFM for corrosion processes monitoring is presented.

**Keywords:** corrosion, atomic force microscopy, synchronous detection, nonlinearity.

### System pomiarowy do nieliniowej spektroskopii powierzchni mikroskopem sił atomowych podczas monitorowania procesów korozji

#### Streszczenie

Oprócz tradycyjnego zastosowania, mikroskop sił atomowych (ang. *Atomic Force Microscope*, AFM) [1] (rys. 1) umożliwia wiele dodatkowych pomiarów, wśród których wyróżnić można obrazowanie powierzchni próbki wyższymi harmonicznymi [6-13]. W trybie półkontaktowym [2-4] w odpowiedzi mikrobekki zauważane są wyższe harmoniczne będące wynikiem oddziaływań nieliniowych [2-6]. Wykorzystanie wzmacniacza fazoczułego (ang.: lock-in amplifier) umożliwia obrazowanie topografii badanej próbki [12, 13] nawet do dwudziestej harmonicznej. Otrzymywane obrazy umożliwiają wyostrenie elementów niewidocznych podczas tradycyjnego skanowania. Ma to związek z wyostroną czułością wyższych harmonicznych na struktury topograficzne. Takie obrazowanie jest szczególnie przydatne do badań próbek niejednorodnych [12], gdy łatwo wyróżnić przez obserwację poziom niektórych harmonicznych [6]. Oprócz wyraźnych zmian przy przejściach między różnymi substancjami, pomiary w cieczach wykazują się zwiększoną wrażliwością na lokalne różnice w elastyczności i interakcji geometrii [10]. Przedstawiono system do obrazowania powierzchni za pomocą wyższych harmonicznych w pomiarach ECAFM (ang. *Electrochemical Atomic Force Microscopy*, ECAFM, rys. 2). Prezentowane rozwiązanie (rys. 3, rys. 4) umożliwia jednoczesną rejestrację obrazów do szóstej harmonicznej podczas monitorowania procesów korozji

(rys. 5, 6). Proponowany system pozwala na uzyskiwanie dodatkowych informacji o zjawiskach korozji zachodzących na skanowanych powierzchniach. Taka metoda stanowi uzupełnienie istniejących metod analizy powierzchni korozyjnych [14-16].

**Słowa kluczowe:** korozja, mikroskop sił atomowych, detekcja synchroniczna, nieliniowości.

### 1. Introduction

An atomic force microscope (AFM) [1] is one of the surface topography imaging methods in the nanoscale. The AFM uses a cantilever beam to place a small scanning tip in close proximity to the sample. The cantilever is a tip-sample force-sensing device [2] which creates a topography map of the specimen. Typically, the force is measured as deflection of the cantilever detected by a laser beam (Fig. 1). This measurement method enhances our ability to characterize materials, map variations in chemical composition, and investigate the properties of nanostructures.

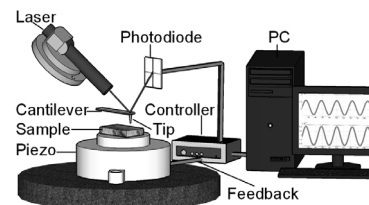


Fig. 1. Schematic illustration of the AFM  
Rys. 1. Schemat budowy AFM

In the standard surface characterization the dynamic AFM (DAFM) is mostly used. Dynamic modes are observed in tapping-mode (TMAFM) where the tip briefly interacts with a sample once per the oscillation cycle. There are two fundamental measurement methods based on TMAFM. The first one is the determination of surface properties by measuring a phase lag between excitation and response signals. However, interpretation of this phase imaging is still discussed. The other one is the detection of higher

harmonics of excitation frequency, present in the response signal [3-5]. Intensities of the higher harmonics correspond to nonlinearity of the tip-sample forces [6].

## 2. Theoretical background

The vibrating cantilever can be considered as a linear system because of its small amplitude of oscillations when compared with its dimensions [7]. The equation of motion of a cantilever tip can be approximated by a one-dimensional forced harmonic oscillator with damping, yielding the nonlinear, second order differential equation [3]:

$$m\ddot{x} + \frac{m\omega_0}{Q}\dot{x} + kx = F_0 \cos(\omega t) + F_{ts}(d_{ts}) \quad (1)$$

where:

- $m$  - mass of the free cantilever,
- $Q$  - quality factor of the free cantilever,
- $k$  - force constant of the free cantilever,
- $\omega_0$  - angular resonance frequency of the free cantilever,
- $\omega$  - angular frequency of the driving force at time  $t$ ,
- $F_0$  - amplitude of the driving force at time  $t$ ,
- $F_{ts}(d_{ts})$  - tip surface interaction force,
- $d_{ts}$  - tip-sample distance.

$F_{ts}$  is a sum of the attractive van der Waals force and the repulsive force due to elastic interactions between the tip and the surface, represented by Lennard-Jones force function. In long-range intermolecular distance  $a_0$ , close to the surface ( $a_0 \leq d_{ts}$ ), the electrostatic force can be neglected. Thus,  $F_{ts}$  is represented only by the van der Waals force. When the tip has contact with the surface ( $a_0 \geq d_{ts}$ ), the elastic interaction of the tip should be included. So that systems without energy dissipation can be modeled using Derjaguin-Müller-Toporov (DMT) [7-10]:

$$F_{ts}(d_{ts}) = \begin{cases} -\frac{HR}{6d_{ts}^2} & \text{for } d_{ts} \geq a_0 \\ -\frac{HR}{6d_{ts}^2} + \frac{4}{3}E\sqrt{R}(a_0 - d_{ts})^{3/2} & \text{for } d_{ts} \leq a_0 \end{cases} \quad (2)$$

where:

- $H$  - Hamaker constant,
- $R$  - the tip radius,
- $E$  - combined Young modulus of the tip (effective contact stiffness) [3]:

$$E = \left( \frac{1-\nu_t^2}{E_t} + \frac{1-\nu_s^2}{E_s} \right)^{-1} \quad (3)$$

where:

- $E_t$  - elastic moduli of the tip,
- $E_s$  - elastic moduli of the sample,
- $\nu_t$  - Poisson ratio of the tip,
- $\nu_s$  - Poisson ratio of the sample.

The Eq.(1) have to be expanded by the additional factor representing a force acting on the cantilever when oscillating at a separation  $d_{ts}$  to the surface in liquid [10]:

$$m\ddot{x} + \frac{m\omega_0}{Q}\dot{x} + kx = F_0 \cos(\omega t) + F_{ts}(d_{ts}) + F_{hyd}(d_{ts}) \quad (4)$$

The nonlinearity of Eq.(2) causes generation of higher harmonics in the cantilever oscillations. The shape of  $F_{ts}$  force function is determined by local physical properties of the specimen. This fact should be helpful in recognition of transition between different parts of the investigated surfaces. Identification of such transitions is especially important during corrosion

process monitoring. Those measurements are performed by electrochemical AFM (ECAFM). That method required a dedicated adapter enabling insertion of the sample in the electrolyte under the controlled potential (Fig. 2).

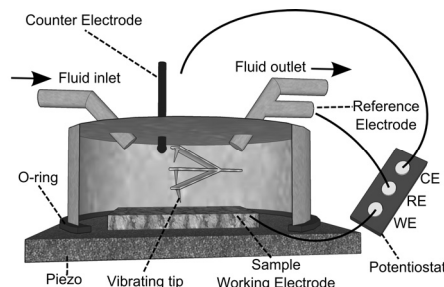


Fig. 2. Schematic illustration of ECAFM

Rys. 2. Schemat budowy ECAFM

## 3. Measurement system

In order to provide higher harmonics measurements, the corresponding signals should be observed. The excitation signal from a generator (GEN) and the cantilever response from a photodiode (DFL) (Fig. 3) are recorded by the data acquisition board. Moreover, it is necessary to record the synchronization signal, which is significant to recognize when every single pixel is scanned. After recording of these three signals for the scanned structure, the data are processed by a dedicated software. During processing the fragments of excitation and response signals are separated for every single pixel. In the next step a software lock-in amplifier to detect higher harmonics was used (Fig. 4). The final image is created by combining the lock-in amplifier calculations.

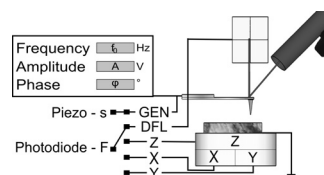


Fig. 3. Parameters of the excitation GEN and measuring response signal DFL

Rys. 3. Parametry sygnału pobudzenia GEN i mierzona odpowiedź dźwigni DFL

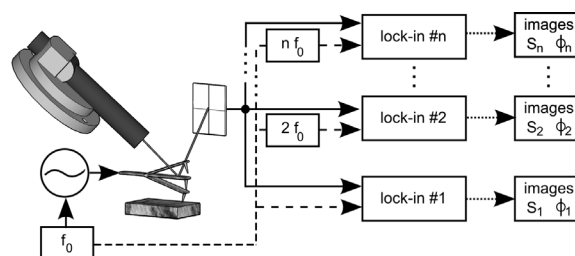


Fig. 4. Higher harmonics detecting by a lock-in amplifier

Rys. 4. Wykrywanie wyższych harmonicznych za pomocą wzmacniacza fazoczułego

The system was described in detail in the previous works [11, 13]. Such measurements require specialized microscopic probes [12] or are limited to an image of the selected harmonic. The presented solution is more general and determines images for a few different harmonics. The previous works [11, 13] proved that the described measurement system worked properly. Fig. 5 shows the sample of glass fiber embedded in resin. Fig. 6 presents the same specimen visualized by higher harmonics images with the area invisible in traditional pictures and marked by the striped line.

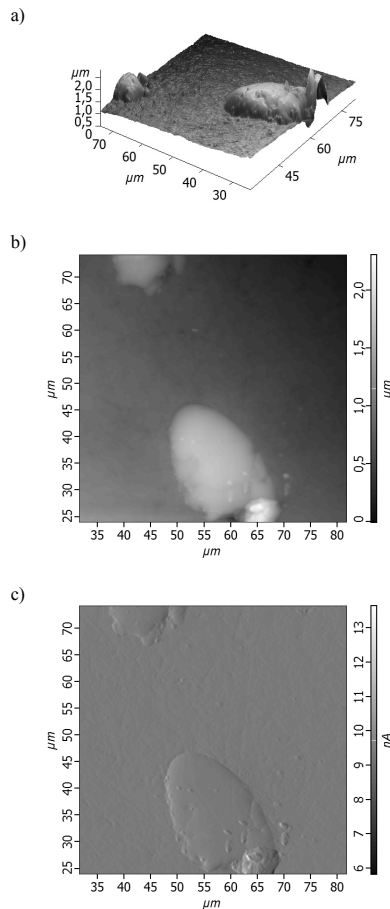


Fig. 5. The sample of glass fiber embedded in resin: a) three dimensional picture, b) topography, c) control error

Rys. 5. Obraz próbki włókna szklanego w żywicy: a) obraz trójwymiarowy, b) topografia, c) sygnał sprzężenia zwrotnego

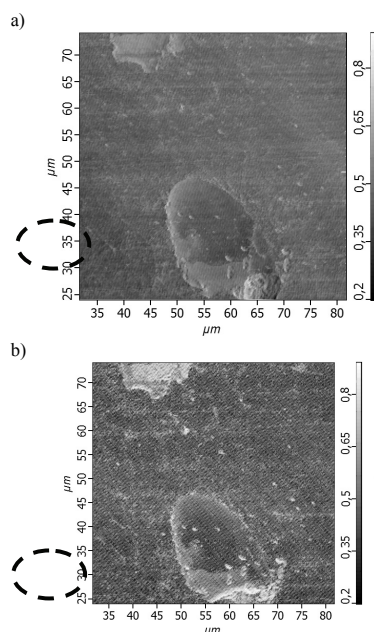


Fig. 6. Higher harmonics images: a) amplitude of the second harmonic b) amplitude of the sixth harmonic

Rys. 6. Obrazy wykonane za pomocą wyższych harmoniczných: a) drugiej, b) szóstej

The proposed measurement technique gives some additional information about physical properties of the scanned surface but there is no clear relation between higher harmonics components and physical properties of the studied structure. We hope that this

method can be successfully applied to characterize a corroding sample but further and more thorough investigation is necessary to establish relations between higher harmonics and physical properties of the investigated sample. Fortunately, we can expect identification by images of higher harmonics of the locally charged area which are caused by varied corrosion intensity and its stage. The charged area will influence the force between the specimen and the tip and therefore will influence the intensity of nonlinear effects. As a result such method would characterize the corroding sample in a more efficient way, which is important in any scientific field related to anticorrosion protection [14-16].

#### 4. Conclusions

The presented results lead to increased functionality of the atomic force microscope by using the developed measurement system. The system can be applied to various atomic force microscope systems at relatively modest costs. Such measurements should provide additional information about corrosion processes taking place on the scanned surface and assure better characterization of their physical properties.

#### 5. References

- [1] Binning G., Quate C. F., Gerber C.: Atomic force microscope. *Physical Review Letters* 56, 1986.
- [2] Garcia R., San Paulo A.: Attractive and repulsive tip-sample interaction regimes in tapping-mode atomic force microscopy. *Physical Review B* 60, 1999.
- [3] Garcia R., Pérez R.: Dynamic atomic force microscopy methods. *Surface Science Reports* 47, 2002.
- [4] Garcia R., San Paulo A.: Dynamics of a vibrating tip near or in intermittent contact with a surface. *Physical Review B* 61, 2000.
- [5] Lee S. I., Howell S. W., Raman A., Reifengerger R.: Nonlinear dynamics of microcantilevers in tapping mode atomic force microscopy: A comparison between theory and experiment. *Physical Review B* 66, 2002.
- [6] Sahin O., Quate C. F., Solgaard O.: Resonant harmonic response in tapping-mode atomic force microscopy. *Physical Review B* 69, 2004.
- [7] Stark R. W., Heckl W. M.: Higher harmonics imaging in tapping-mode atomic-force microscopy. *Review of Scientific Instruments* 74, 2003.
- [8] Stark R. W.: Spectroscopy of higher harmonics in dynamic atomic force microscopy. *Nanotechnology* 15, 2004.
- [9] Hillenbrand R., Stark M., Guckenberger R.: Higher-harmonics generation in tapping-mode atomic-force microscopy: Insights into the tip-sample interaction. *Applied Physics Letters* 76, 2000.
- [10] Preiner J., Tang J., Pastushenko V., Hinterdorfer P.: Higher Harmonic Atomic Force Microscopy: Imaging of Biological Membranes in Liquid. *Physical Review Letters* 99, 2007.
- [11] Babicz S.: The study of harmonic imaging by AFM. *Pomiary Automatyka Kontrola* 12, 2011.
- [12] Sikora A., Bednarz Ł.: The implementation and the performance analysis of the multi-channel software-based lock-in amplifier for the stiffness mapping with atomic force microscope (AFM). *Bulletin of the Polish Academy of Sciences: Technical Sciences* 60, 2012.
- [13] Babicz S., Zieliński A., Smulko J.: Enhancing capabilities of Atomic Force Microscopy by tip motion harmonics analysis. *Bulletin of the Polish Academy of Sciences: Technical Sciences* 61, 2013 (during the issue).
- [14] Hoja J., Lentka G.: Analyzer for spectroscopy of low impedance. *Metrology and Measurement Systems XVI* (1), 2009.
- [15] Niedostatkiewicz M., Zielonko R.: Investigation on accelerated impedance spectrum measurement method with mulisine signal stimulation objects. *Metrology and Measurement Systems XVI* (4), 2009.
- [16] Antoniuk P., Strąkowski M., Pluciński J., Kosmowski B.: Non-Destructive Inspection of Anti-corrosion protective coatings using optical coherent tomography. *Metrology and Measurement Systems XIX* (2), 2012.